

# Datasheet




Mouse mAb to **CD10**  
Clone **FR4D11**  
Isotype **IgG1-κ**

## Source

A BALB/c mice was immunized with Raji cells.  
Fusion partner: X63-Ag8.653.

## Specifications

FR4D11 reacts with high affinity to CD10 or CALLA, a cell surface enzyme with neutral metalloendopeptidase activity, inactivating a variety of biologically active peptides. CD10 is a 100 kDa glycoprotein, expressed on 70% of pre-B ALL-cells (common ALL), but also on early lymphoid progenitor-cells in bone marrow and fetal liver. Other normal CD10 positive tissues include renal epithelium, fibroblasts and germinal centre B-cells. Density of CD10 antigen has been shown to be related to cell differentiation and may have prognostic value for B-cell lineage acute leukemia. CD10 is also present on breast myoepithelial cells, bile canaliculi, fibroblasts, with especially high expression on the brush border of kidney and gut epithelial cells.

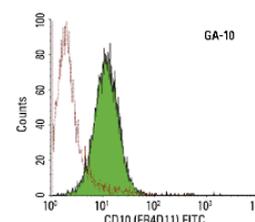

**Figure 1:** Burkitt's lymphoma stained with FR4D11 (FACS).

## Species reactivity

Positive: human.

## Applications

FR4D11 can be used for the classification of acute leukemias and childhood ALL prognosis (patients (CD10+ have a better prognosis than CD10-).

Flow cytometry	Frozen sections	Immunofluorescence
+	+	+

## Format

Produced in tissue culture, contains no host Ig. Antibodies are affinity purified and presented in PBS with 0,02% sodium azide.

Stored at 4°C-8°C, shelf life is at least 24 months after purchase.

## Dilution advice

- Flow cytometry (0,5-1,0 µg/million cells in 0,1 ml).
- Immunofluorescence (0,5-1,0 µg/ml).
- Immunohistology (1-2 µg/ml for 30 min at RT; an appropriate antigen retrieval method for staining of formalin-fixed tissues has not been established to date).

## Positive control

Raji cells, tonsil, small intestine or kidney.

# Datasheet



## References

- Brown, B. et al., *J. Natl. Canc. Inst.*, **55**: 1281-1289 (1975).
- Tran-Paterson, R. et al., *Blood*, **76**: 775-782 (1990).
- Doerken, B. et al., in Knapp, W. et. al. (eds), *Leucocyte Typing IV*, Oxford Univ. Press, pp 33-34.
- Lavabre-Bertrand, T., et. al., *Cytometry*, **18**: 209-217 (1994).